

MADE IN INDIA



UNIK SURGICAL PVT. LTD.

Manufacturer, Importer, Exporter & Suppliers

Medical Furniture, Lab Furniture, Non Medical
Furniture, MOT, MGPS, Laundry Equipments, Kitchen Equipments, CSSD, TSSU,
BMW, Medical Equipments, Mortuary Equipments, Nurse Call System,
Lab Equipments & General Orders Supplier



MEDICAL GAS PIPELINE



QUALITY | TECHNOLOGY
COMMITMENT | SERVICES



UNIK SURGICAL PVT. LTD.

Founded in 1999, Unik Surgical Pvt. Ltd. has been at the forefront of healthcare infrastructure development in Eastern India. What began as a small enterprise has grown into a leading supplier and manufacturer of medical and surgical equipment and furniture.

With our headquarters in Patna, Bihar, a regional office in Greater Noida, Uttar Pradesh, and a manufacturing facility in Chapra, Bihar, we serve clients across the region with unmatched expertise and a strong commitment to quality. Our journey is defined by continuous innovation, strict adherence to international standards, and an unwavering dedication to improving healthcare facilities across India.

Contracting

Project Report

Layout

Configuration

Consulting

Installation

On-site Training

Startup

Costing

Technical Support

ABOUT DIRECTOR

Manoj Kumar Singh is a healthcare industry leader known for his practical approach and deep commitment to improving patient care. With years of hands-on experience in medical equipment manufacturing, he founded Unik Surgical with a simple goal — to make reliable, high-quality surgical solutions accessible to healthcare facilities of all sizes.

What began as a trader in Patna has now grown into a trusted name, serving hospitals and clinics across the country. Manoj's focus on precision engineering, transparent business ethics, and long-term client relationships has helped Unik Surgical earn the confidence of medical professionals.

His leadership continues to guide the company toward innovation and responsible growth, while staying rooted in the belief that every piece of equipment should ultimately contribute to better patient outcomes.



International & National Accreditation of Bodies



एन. एस. आई. सी.
NSIC
ISO 9001 : 2008

CERTIFICATIONS



QUALITY | TECHNOLOGY
COMMITMENT | SERVICES

MEDICAL GAS OUTLET POINTS

BS Standard (Model: Amber)

Front Loaded with Probe



Features

- Complies with BS EN 737-1, BS 5682 and HTM 2022 standard regulations.
- CE approved.
- The Medical Gas outlet combines ease of installation and low cost maintenance.
- Available for head wall unit/wall/pendant units.
- Double locking (secondary) facility.
- Front fascia of Gas outlet are non-interchangeable and gas colour coded.
- Secondary isolating valve allows the maintenance of the gas outlet without the need for a pipeline system shutdown.

British Standard (BS) Medical Gas Outlets are fixed to the gas source by safety nut preventing non-competent personnel to connect the outlet to the wrong gas source.

The releasing rings are gas specific via color and different shapes and it is absolutely impossible to mix them. Geometrical and dimensional connector keying according to gas used safety. Non-return valve for no interruption in gas supply. Plastic fool-proofing design to prevent gas interchangeability during reassembly. Color coding for different gas source and only accepts British compatible adapters. It can be used on O₂, Compressed Air, Vacuum and N₂O.

PB Standard (Model: Sapphire)

Front Loaded with Probe



Complies with NFPA-99

Features

- Flow rate exceeds the requirements of NFPA.
- Accept only puritan-Bennett gas adapters.
- 100% Pressure testing of each outlet and cleaned for medical gas service.
- Indexed to eliminate interchangeability of gas services.
- Up to 1" [25mm] wall thickness adjustment.
- 360° Swivel inlet Pipe for easy installation.

A large color coded front plate shall be used for ease of gas identification and aesthetic appeal.

It is manufactured with a 6-1/2" [165 mm] length type "K" 1/2" [12.7mm] outlet diameter (3/8" nominal) size copper inlet pipe stub, which is silver brazed to the outlet body. The body shall be of 1-5/16" [33mm] diameter, one piece brass construction. For positive pressure gas services, the outlet shall be equipped with a primary and secondary check valve and the secondary check valve shall be rated at a

maximum of 200 psi [1,379 kpa]. In the event the primary check valve is removed for maintenance, outlets body shall be gas specific by indexing each gas service to a gas specific dual pin indexing arrangement on the receptive identification module.

Model: Parkodex with Probe

(Hook Loose Type)



Features

- The Medical gas outlet combines ease of installation and low cost maintenance.
- Non-interchangeable self-sealing quick connect gas outlet for O₂, N₂, Air, CO₂.
- Available for head wall unit/wall/pendant units.
- Double locking (Secondary) facility.
- Front fascia of gas outlet are interchangeable and gas colour coded.
- Secondary isolating valve allows maintenance of the gas outlet without need for a pipeline system shutdown.

DIN Standard (Model: Jasper) with Probe



Features

- Geometrical and dimensional connector keying according to gas used.
- Two-part construction for simple and fast separation from base block.
- Double locking prevents loss of gas due to leakage and for facilitating online servicing also.
- Gas-specific block and socket prevent any risk of false operation as interchanging of gas service and non-competent personnel mix connecting the outlet to the wrong gas source.
- Equipped with non-return valve for the convenient installation and servicing replacement possible without blocking the station.
- Push to release mechanism adaptor will automatically click back into first click lock position. Gas flow to the equipment is interrupted by pressing release ring on the fascia.
- Available for different kinds of installation, e.g. in walls, consoles, bed head trunkings, ceilings or pendants.
- Module design minimize the number of joints where dust and bacteria can collect. Thus it is very easy to clean and maintain.
- Worldwide color-coding and lettering identification for clearly indicating the specific gas.

German Standard Medical Gas Outlets are designed for quick connection with no risk of gas cross connection for gas distribution in centralized areas of hospitals, dental clinics and veterinary facilities; patient rooms, operating areas, surgery centers and dental offices.

OXYGEN FLOW METER with HUMIDIFIER BOTTLE

BPC Flow Meter



Features

- It will be made of nickel - plated brass body.
- It will be backpressure compensated, double jacket with flow gradations marked on it. Desired flow will be settled through rotatory knob.
- Flow range should be between 0-15lpm.
- Transparent and unbreakable type humidifier bottles, to provide bubbletype humidification, made of polycarbonate to be installed with nut and nipple type connections.
- Humidifier bottle is available with safety valve features.

Oxygen flow meter is a pressure compensated flow meter that delivers an accurate prescribed flow to patients from the gas source. Our oxygen flow meter provides an accuracy of $\pm 5\%$ of full scale reading and is pressure compensated to ensure that the flow meter maintains its accuracy independent of varying gas supply pressures. It is suitable for most type of respiratory therapy.

THEATRE SUCTION TROLLEY

Theatre Vacuum Unit Modle JK /TVU/01



Theatre Vacuum Unit Modle JK /TVU/02



It will be trolley mounted. The unit will include one regulator having gauge mounted on the stand. Having two reusable each 2000ml plastic collection bottles mounted on base of trolley unit connected with regulator and low pressure tube inter connections. The base of the trolley is fitted on castor wheel arrangement to move the system from one place to another place easily.

VACUUM REGULATOR



Features

- Meets EN 739 standards.
- Provides an adjustable, continuous vacuum level up to 760mmHg.
- REG/OFF switch allows the quick restoration of a pre-adjusted vacuum level.
- Continuous suction controller guaranteeing reliability and fine adjustment.
- Plastic hosing and cover preventing from risk of corrosion.
- REG/OFF indicator on the side allows the care giver to operate the vacuum regulator easily.
- Non Interchangeable Screw Threads (NIST) guarantee easy and safe operation.
- Precise vacuum control with non-slip hand wheel.
- Lightweight plastic construction
- Colour coded body for additional safety
- User-friendly gauge marked for high, medium and low pressure.
- Multiple connections and mounting options.
- Available with overflow safety trap.

2-stages Continuous Model Vacuum Regulator is to measure and adjust the vacuum level within the context of medical suction. It is used to remove substances from the patient's body during surgical and other procedures. The vacuum regulator should be connected to a wall source of vacuum using a direct probe or through a rail mounting system. Applications include all medical and emergency applications.

VACUUM JAR



Features

- Clear polycarbonate is non-toxic, exhibits low protein binding and is easy to clean.
- Screw-fit plastic cap and float assembly can be steamed or sterilized in autoclave.
- Cap fitted with float safety valve automatically and efficiently shuts Vacuum fully when it is full.
- Cap inlet and outlet tubing connectors molded-in unique graphic indication of patient/vacuum for easy-use.
- Molded ergonomic handle enables it to hang it directly onto the wall side which reduces the risk of fluid contact and makes transporting simple and efficient.
- O Ring over the top of bottle for a positive, leak proof seal without generating any noises or vibrations harmful to the patient.
- Available in 600ml, 1000ml and 2000ml.


Feature in 2000ml Jar

- In 2000ml capacity, the cap/lid is clipped on to avoid the fluid spread outside the bottle at the time of opening the Jar lid. This helps in easy removal of the lid from Jar.

Reusable Collection Bottle is used to collect liquids or mucus resulting from medical suction. Constructed of polycarbonate, it can be sterilized in autoclave countless times without losing quality even at 134° C for 15mins. It can withstand hard usage and provides trouble free service over a long time.

AVSU ZONAL VALVE BOX

Zone Valves and Zone Valve Boxes provide versatility and value through medical gas or vacuum pipeline flow prevention.

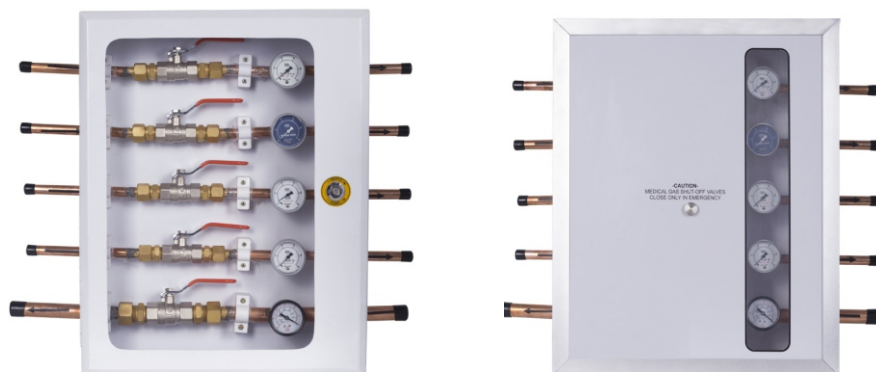

Features

- Offered in single (houses 1 valve) or multiple (houses 2 to 6 valves) types.
- Features a pull-out, removable opaque window.
- Adjustable to different wall thicknesses.
- Include color-coded gas identification labels.
- Blowout proof stem.
- Available in locking or non-locking.
- Swing-out design to prevent damage during brazing operation.
- 100% hydrostatically tested.
- 3-piece ball type design for ease of maintenance.

MANIFOLD SYSTEM (OXYGEN/N₂O)

The manifolds manufactured by us are extendable & Non-extendable type with choice of size of cylinders bank for bulk 'D' type jumbo cylinders.

Manifold is suitable to withstand a working pressure of 140Kg/cm³.

Manifolds have high pressure copper annealed tail pipes with one end having Brass adaptor suitable for cylinders. Manifolds have Brass Non Return Valves (NRV) for individual cylinders (Cylinders to be provided by the client).



AUTOMATIC GAS CONTROL PANELS (DIGITAL/ANALOG)

Switchover gas manifolds provide a standby gas supply from two equal banks of cylinders, with automatic changeover. The manifold is convenient through regulating gas cylinders pressure to provide uninterrupted gas. It provides a cost effective solution to trailing hoses around the operating room.

The automatic changeover is done through a Fully Automatic Gas Control Panel.

Spectrolite



Features

- Complies with NFPA 99.
- Fully automatic with manual reset.
- Lockable metal outer cover provides durability and fire-resistance enhanced safety.
- Automatic switch by saddle valves between two gas sources and is used for refilling gas.
- Pressure-reducing design maintains a reliable stable output pressure.
- Pressure gauge/LED display shows clearly the pressure at each gas source and the total pressure at the outlet.
- Full bore valves for minimum pressure loss can be locked in open or closed positions.
- Provides maintenance or major works through isolation facilities of gas distribution system.
- Microprocessor based PLC for ease of maintenance and error reporting.
- Gas specific color-coded and non-interchangeable hoses connections results in enhanced safety.
- Audio and light alarms turn on automatically when one of the gas sources is exhausted.
- LED display readable even in poor lighting connections.
- Filter screens out all the impurities in the gas sources.
- User-definable hose length upon requests.
- Remove status panel indicates local status.
- Manifolds available for medium and high flow applications.
- Passes 100% pressure tests.
- Available in handling gases of O₂, Air, N₂O and CO₂. Input power 110V AC to 240V AC, 50 to 60Hz.

SEMI AUTOMATIC

We manufacture semi automatic gas control panels with manifolds. These equipment are used for checking the gas level. Fitted with control panels, our product are known for their ability to yield accurate results.




POWER COLUMN

Power Column samples the organization of medical service equipment in one aesthetically appealing vertical headwall.

It is a modular designed product to provide a clean, sanitary and ordered environment. Furthermore it can enhance the quality in care and medical environment. Due to its systematized design, a network provided for medical care service is applicable to the use of secondary instruments and assemblies.

The main structure consists of fluorescent panel, aluminum component and side track. Also includes power lines and medical gas piping. Side tracks are installed at four at angles, and the length is equal to the column. All the assemblies required by medical care personnel can be easily and securely placed on the side track of this column. Applications in mid to acute level environment.

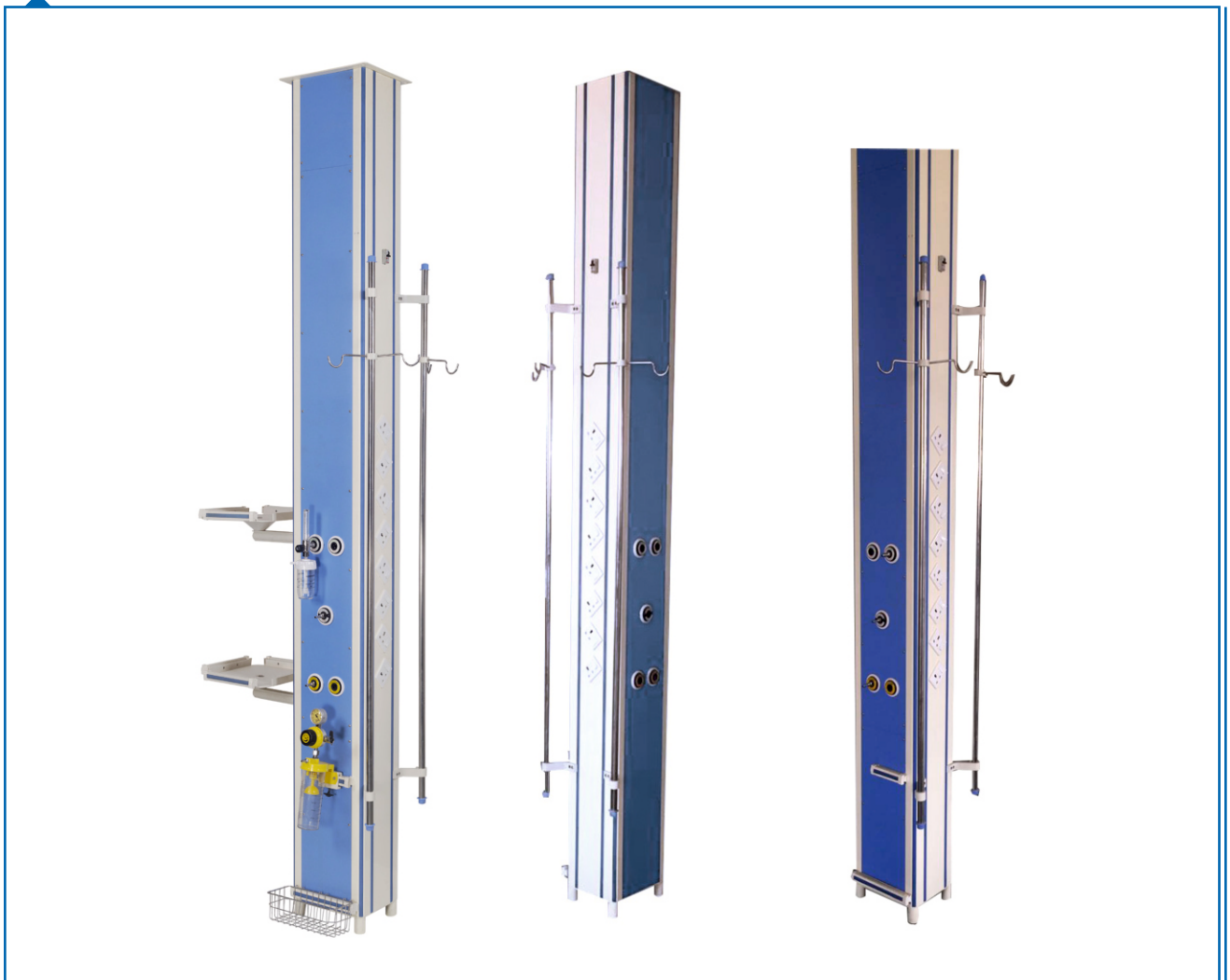

Features

- Designed for ease of maintenance and infection control.
- Built as a four-sided unit that offers 360° access to services and equipment management for critical care areas.
- Can be customized and upgraded to meet your facility's needs.
- Offered in a variety of widths and heights, tailored to meet your space needs.
- Wide selection of laminate colors are available.
- Multiple rail system allows for superior accessory management.
- Removable laminate panels allow for simplified maintenance.

BED HEAD PANEL

Bed Head Units/Bed Head Panels are constructed from aluminum extrusion profiles with round and smooth shape corners and shall be in different sizes & length per each patient. The unit shall have three physically separated compartments (for low and extra low voltage, electrical services and for medical gas services). The unit shall provide a safe, efficient means of delivering services to patients/staff in both general and special care applications.

The unit shall have smooth curved surfaces, no visible screws. Installation of the unit shall be carried out easily by aluminium rail-brackets which allow simple horizontal and vertical alignment taking into account any wall high spots or unevenness. Removal of the covers for maintenance and pipeline connection shall be easy.

The unit shall consist of a number of modular extrusions capable of being joined together to form a carcass to suit the particular application. The unit shall have a backspin which shall incorporate an area which is for medical gas pipelines.

Segregation of services ie. Extra low voltage (<48 volts), low voltage (<500 volts) and medical gases shall be maintained throughout. Medical Gas and Electrical services shall be separated by inside unit. Electric power feedings shall be 220V/50Hz. The middle compartment of the unit shall be used for medical gas outlets and electrical services but gas services and electric services shall be physically separated. Upper and lower compartments of the unit shall be used for reading through indirect lights and necessary electrical wirings. The lower cover of the unit shall have compartment for reading through light illumination. The upper cover of the unit shall have two openings for indirect light illumination through a transparent sheet. Both sides of the unit shall be closed by aluminum covers. All cables and terminal blocks used inside the unit shall be non-flammable type.

The unit shall be manufactured in accordance with international and European standards, ISO 11197 and EN 793 for "Medical Supply Units" to HTM 2015 and 2022 and to latest edition of IEE and NEC regulations. The aluminum extrusion unit shall have an electrostatic powder paint finish on all external surfaces. Unit shall be supplied pre-piped, pre-wired and fully tested. The minimum mounting height of the unit should ideally be 48 inches (4') from the finished floor level. All removable covers of the unit shall be earthed. Electric power feedings shall be 220V/50Hz.

Lumipan H



The standard range of Patient Bed Head Unit shall be "CE" marked under the Medical Devices Directive 93/42/EEC. Under this directive, the specified products shall be classified as Class IIa Medical Devices. (We have different types of models and specifications can vary for each model)


THE UNIT SHALL CONSIST OF FOLLOWING PROVISIONS.

1. Medical Gas outlets	5. Infusion pole	9. Reading lamps
2. Electrical Switch Sockets	6. Vacuum jar holder	10. Medical rails
3. Low voltages	7. Monitor Trays	
4. IV Hook	8. Up and down lights	

Silk H



Critrunk H

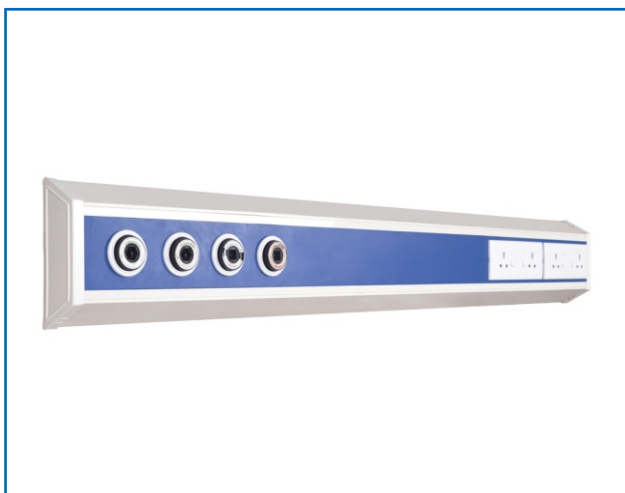

Lumipan D



Robustrunk



Starlite H



Silk Smart V



Spinel V



Lumipan V



Citrunk V



Silk Smart V



OPERATION THEATRE SURGICAL PENDANT

The Ceiling Pendant Systems shall provide convenient positioning of medical equipments, medical gas outlet points, electricals and speciality services. The pendant shall be mounted to the ceiling with extremely robust carcass system comprised of stainless steel angles and supporting rods. The Ceiling Pendant shall comply with HTM2022. The support arms shall be extremely robust and revolve on high quality bearings, so that the pendant head glides smoothly and quickly to any desired position. The multipurpose pendant can have single arm or double arm. The arms shall be available as follows:

-1000mm movable arm with 600mm of vertical movement arm shall be capable of 330 degrees of rotation and be fitted with stops, which can be easily adjusted to suit the desired mode of operation. The arm may be fitted with pneumatic brakes to prevent inadvertent movement. The pendant service heads shall be modular. The head shall be capable of accepting a range of shelves, accessories rails or other accessories.

The Pendant Service Heads shall be in 500mm length and supplied with medical gas outlets, BS 13 amp unswitched sockets. Medical gas outlets and probes shall fully complying with BS 5682:1992. Gas connection between pendant and pipeline shall be carried out via antistatic medical gas hoses comply with EN 739:1998. Both ends of the hoses shall be connected by gas specific NIST (Non-Interchangeable Screw Thread) connectors. The heavy duty electric arm shall be capable of supporting 90 kg load at maximum reach. Pendant arm shall be RAL 9002 electrostatic powder painted. The medical rail poles shall have an anodized finish.

The pendant shall have an electric junction box mounted to the carcass system over false ceiling for connection of incoming electrical services: segregated into three compartments, electrical, communications (low-voltage), and critical branch power, as required. Preinstalled electrical equipment shall be wired to the junction box. The pendant shall have ground bar located in junction box. Ground bar shall be copper with compression screws, which shall accept 4 AWG, or smaller building service ground wire. All cables used inside the pendant shall be non-flammable low smoke halogen free type. Cables shall be passed through flexible spiral conduit. Surgical Pendant shall be equipped with the following services:

1. 2 no of Oxygen (O₂) outlet
2. 1 no of Nitrous Oxide (N₂O) outlet
3. 1 no of Medical Air 4 bar outlet
4. 1 no of Medical Air 7 bar (Surgical Air) outlet
5. 2 no of Medical Vacuum outlet
6. 2 no of electrical box made from moulded white self-extinguishing ABS at both sides of the console
7. 4 no of double unswitched BS 13A standard electrical socket (UPS)
8. 8 no of potential earthing nodes
9. 1 no of brake button for the control of pneumatic brake
10. 2 no of aluminium extrusion accessories rail
11. 1 no of shelves in 500x470x25 mm sizes
12. 4 no of Pressure gauges for oxygen, nitrous oxide, air 4 and air 7 bar
13. 1 no of vacuum gauge

The standard range of Operation Theatre Pendant shall be "CE" marked under the Medical Devices Directive 93/42/EEC with approval body number. Under this directive, the specied products shall be classied as Class IIa Medical Device

S-2 (V) 600



S-2 (V) 800



OPERATION THEATER PENDANT

S-2 (V)



S-3 (V) 800



S-4 (V) 800



S-3 (H)



S-I (V) 1000



MEDIA BRIDGE (OT/ICU BRIDGE)

We bring you the latest technology in the form of a Media Bridge. Media Bridge is a replacement for pendants. It is not only better economically by providing 360 degree access but it also helps in maintaining the class of clean room specified as there is no obstruction in the laminar flow of air.



Features

- Aluminum Extrusions: Robust construction design with all suitable segregation.
- Medical gas outlet points (all the four sides of media bridge).
- Electrical switch sockets 5/15amp (All four sides).
- Low voltage segregation, RJ 45 sockets.
- Manually Horizontally (Left to Right) Movements of monitor trays (Pendants).
- Monitor trays swivel movement up to 330 degree.
- Hanging Arrangements for Infusion pump, IV and other suitable equipments.
- Suitable High pressure antistatic and antibacterial hoses.
- Concealed LED Lights for better luminosity.
- Height adjustable trays (Pendants).

S-4 (V) 600



Movable

S-5 (H)



Laparoscopic
Movable

MOTORIZED PENDANT




Specifications

- Material High strength aluminum alloy or customized
- length of the arms; One: 750mm; Two: 750mm
- Loading capability One 150kg
- Medical gas outlets Oxygen (2 nos.), Vacuum (1 nos.) and Air (1 nos.)
- Equipment shelf 2 pcs or customized
- Equipment drawer 1 pcs or customized
- Infusion sets Include sets sockets, infusion pump, infusion poles. etc
- Power Input 220V; 50-60Hz or customized
- Gas outlets standards: PB/BS/DIN/Parkodex
- Rotary degree 330 degrees in each joint
- Ceiling mounted according to customers' requirements
- Power sockets 6 pcs
- Communication system RJ45 (1 Pcs); RJ15 (1 Pcs)
- Finishing: Antibacterial Powder.

Single Arm Motorized Medical Pendant is mainly used in the Operating Rooms, Intensive Care Units, Emergency Rooms, etc. to supply Medical Gases, Electrical Power, Communicating Ports, and to load the medical equipments such as patient monitors, infusion pumps, surgical instruments, endoscopy equipments, etc.

S-5 (H)

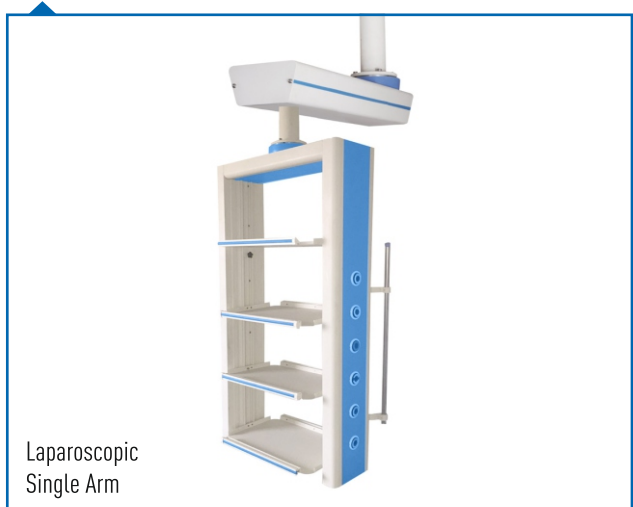

S-4 (V) 800



MEDICAL GAS ALARM (DIGITAL)

The Compact Digital Alarm system shall be complete with a one-year warranty. Each module of the Compact Alarm shall be microprocessor based with individual microprocessors on each display and sensor board. The sensors shall be capable of local (in box) or remote mounting. Each sensor and display unit shall be gas specific, with an error message display for an incorrect connection.

Compact Alarm can accommodate up to 07 gases.

Each specific Area service for the Compact Alarm shall be provided with an LED digital read-out comprising of 0-250psi [0-1,724 kPa] for pressure and 0-30" Hg [-100-0 HgkPa] for vacuum. The digital read-out shall provide a constant indication of each service being measured. It will indicate a green "NORMAL" and a red "HIGH" or "LOW" alarm condition. If an alarm occurs, the "RED" alarm light shall flash and the audible alarm will sound. Pushing the "ALARM MUTE" button will cancel the audible alarm, but the unit will remain in the alarm condition until the problem is rectified.

The default set point on the Compact Alarm shall be $\pm 20\%$ variation from normal condition. In the calibration mode the following parameters shall be field adjustable: High/Low set-points, Imperial/Metric units, Repeat alarm Enable/Disable. Set points shall be adjustable by two on board push buttons.

Each Compact Alarm module shall be microprocessor based and be adjustable.

A repeat alarm function shall, when enabled on the Compact Alarm, be capable of turning on the buzzer again (after a preset time) if the fault condition has not been rectified.

The box shall be fabricated from 18 gauge [1.3 mm] steel with 1/4" [6.4 mm] I.D. type "K" copper pipe for connection to the service line. The box mounting brackets shall be adjustable to accommodate for different wall thickness.

Input power to the alarm is: 115 VAC to 220 VAC, 50 to 60 HZ. Complies with NFPA-99.

All the gases inlets provide the NIST connections.



Features

- Digital sensors can be mounted locally or remotely utilizing.
- True digital LED read-out with red "alarm", normal green trend indicator for each service.
- Psi, kPa or BAR read-out (switch selected).
- High/low alarm set-points are field adjustable for each gas service.
- Repeat alarm, adjustable 1 to 60 minutes.
- Up to 07 gases in a standard configuration.
- Alarm buzzer in excess of 90 decibels.
- Individual microprocessor for each display and sensor module.
- All modules mounted on a hinged frame for easy accessibility.

MASTER ALARM

Master Alarm System is mainly used in the system for providing information and current status of multiple gas alarms and is capable of supporting various gas alarms of different gas categories. The main feature is for remote monitoring of various gas alarm installed in hospital's rooms with displayed color LED's. The main benefits of this product are:

Main Master Alarm

20 Windows Master Alarm

Master Gas Alarm



Features

- Support for multiple gas alarms (Max up to 20).
- Different category gas alarm support (2 Gas to 6 Gas Type).
- Real time and simultaneous monitoring from one single location
- Interactive color LED's push buttons and 7-segment display.
- Buzzer provided for sound alarm.
- Mute and test buttons.

MEDICAL GAS ALARM (ANALOG)




Features

- Pressure sensors can be mounted locally or remotely utilizing.
- True analog gauge read-out with red "alarm", normal green trend indicator for each service.
- True analog gauge read-out with red "alarm" normal green trend indicator for each service.
- High/low alarm set-points are held adjustable for each gas service.
- Repeat alarm, adjustable 1 to 60 minutes.
- Up to 07gases in a standard configuration.
- Alarm buzzer in excess of 90 decibels.
- Individual microprocessor for each display and pressure sensor module.
- All modules mounted on a hinged frame for easy accessibility.

NURSE CALL SYSTEM

Nurse call system acts as an interface between patients and the hospital administration/nurses.

The main unit is placed with the administrator and its subordinate 'Corridor panels' placed outside which are used by the patients to call the nurse. Nurse Call system makes the patient administration efficient and greatly simplifies patient care. The patients can press the appropriate button when he requires assistance, which will then show on the corridor and main panels. If the nurse is near the room, she will see the corridor panel light and immediately assist the patient. Otherwise, the attendant at the main panel can acknowledge the call, which will light up the LED 'waiting' on the bed panel informing the patient his call has been acknowledged. The main panel displays the corridor wise and room wise segregation of panels. The usual problems associated with installing such systems are extensive and complicated wiring. The usage of single cable wiring in Nurse Call System makes sure wiring is easy and not expensive.

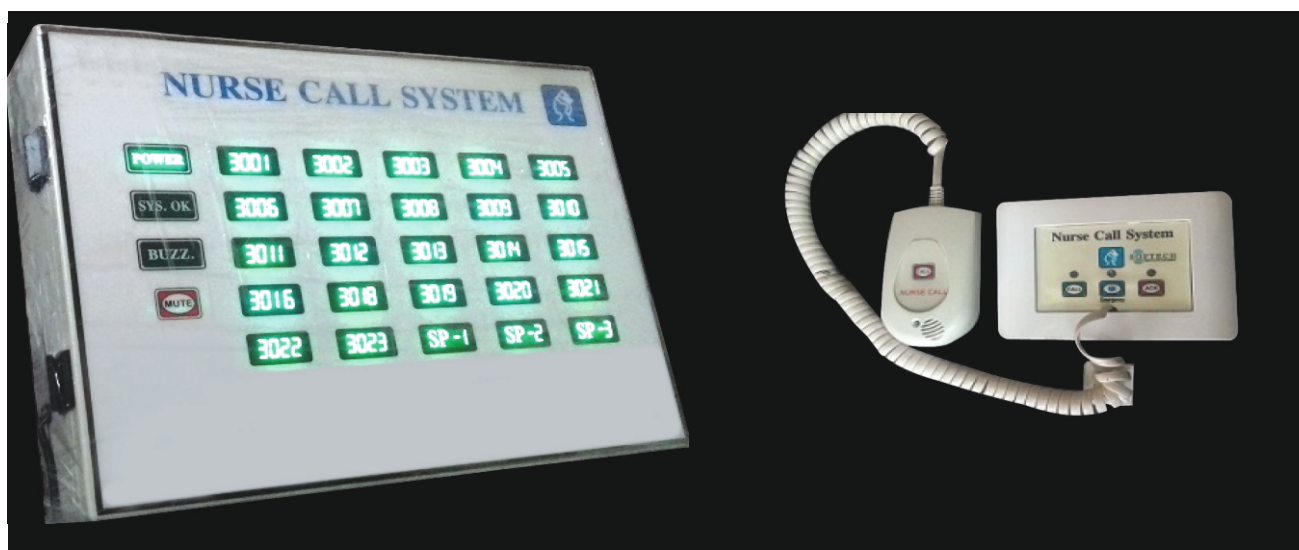



UNIK SURGICAL PVT. LTD.

Head Office : Ground Floor, Angad Apartment,
GD Mishra Path, New Patliputra Colony, Patna - 800013, Bihar (India)

Phone : 0612-3562111

Email : info@uniksurgical.com | uniksurgicalpvtltd@gmail.com

Website: www.uniksurgical.com

Authorized Dealer